ViiSona
ultrasound & education

PORTABLE VETERINARY ULTRASOUND


MULTI-SPECIES FIELD ULTRASOUND

iScan 2 MULTI

As versatile as you are

One rugged, lightweight scanner with interchangeable probes for cattle, horses, swine and small ruminants. A broad 2–14 MHz range, a bright 7″ screen and up to 7 hours per battery make it ready for any call.

- **Interchangeable probes, auto-recognized**
- **Up to 7 hours per battery**
- **Back-fat probe for carcass evaluation**

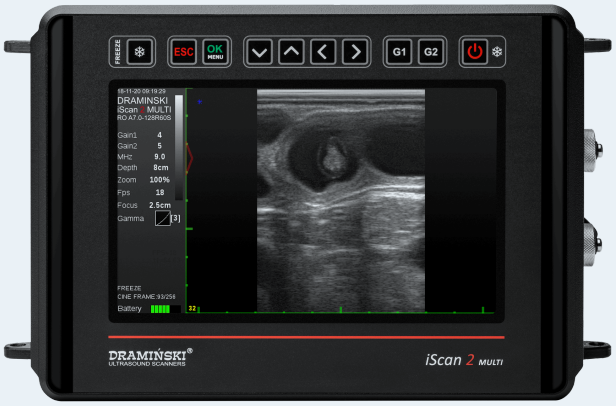

STARTING AT

iScan 2 MULTI scanner
Probes $4,739–$6,495 each

$10,295

2–14 MHz
Frequency range (probe dependent)

7 h
Operating time per battery

2.64 kg
With probe and battery

7″
LCD LED, 500 cd/m²

www.viisona.com · 1-800-672-9182

DRAMIŃSKI TECHNOLOGY

Ideal for

Cattle, horses, goats, sheep and pigs.

ADDITIONAL APPLICATIONS

Mixed large-animal practices that need one scanner for many species.

01	Reproduction and pregnancy diagnosis
02	Fetal aging (cows, horses, sheep, llamas, sows)
03	Back-fat and carcass merit measurement
04	Abdominal examinations
05	Lung examination
06	Ovum pick-up (OPU)

One scanner, many species

Interchangeable probes

Choose from rectal, abdominal and specialized probes on a single port with automatic probe recognition. Swap probes in seconds between patients.

Large, bright display

A 7″ LCD LED screen, 500 cd/m², with a wide viewing angle that stays readable at the chute.

Fully digital imaging

B, B+B and B+M modes with LuciD image enhancement, zoom from 60 to 200%, negative view and split screen.

Complete measurement set

Distance, area, volume, grid and back-fat tools, plus fetal aging tables for cows, horses, sheep, llamas and sows.


Rugged and easy to clean

A duralumin body, splash-proof keyboard and a key layout that suits both left- and right-handed users.

Long working days

Up to 7 hours per battery with up to a 4-hour charge. Store 200 images and 200 cine loops; export to a USB drive or send images to the mobile app with the Wi-Fi antenna.

DIAGNOSTIC IMAGES

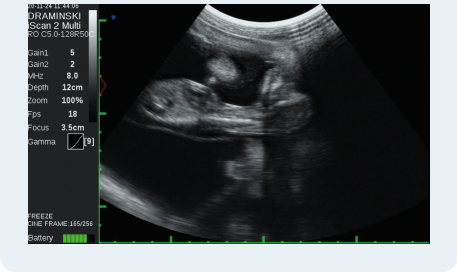

Fetal image on iScan 2 MULTI



Rectal, abdominal and specialized probes



Swappable external battery

Build your configuration

Prices start at the amounts below. Your final price depends on the probes and accessories you choose. Ask us for a quote built around your practice.

ITEM		STARTING PRICE
iScan 2 MULTI scanner (probe sold separately)	Choose one or more probes	**$10,295**
Interchangeable probes	Per probe, depending on model	**$4,739–$6,495**
Wi-Fi antenna, stand and accessories	Priced to your configuration	**Quoted**

Probe options

17 compatible electronic probes, 128 crystals each, recognized automatically. Select at least one probe with your scanner.

RECTAL

163 · Rectal linear, 60 mm face — 4–9 MHz
Cattle and horse repro, pregnancy, fetal sexing; lungs, udder

212 · Rectal convex, R60 mm — 2–8 MHz
Repro, pregnancy, lungs, udder; small ruminants; abdominal organs

191 · Rectal linear (sheep and goats), 40 mm face — 5–10 MHz
Rectal exams in sheep and goats

192 · Rectal microconvex (P-probe), R20 mm — 4–9 MHz
Rectal exam of the horse musculoskeletal system

ABDOMINAL CONVEX

172 · Abdominal convex, R60 mm — 1–6 MHz
Deep abdomen; small ruminant pregnancy; horse musculoskeletal

161 · Abdominal convex, R50 mm — 2–5 MHz
Small ruminant pregnancy; abdomen; horse musculoskeletal

182 · Abdominal convex, R60 mm — 2–5 MHz
Small ruminant pregnancy; abdomen; horse musculoskeletal

156 · Abdominal convex, R50 mm — 2–8 MHz
Small ruminant pregnancy; cattle and horse abdomen and thorax

225 · Abdominal convex, R50 mm — 3–7 MHz
Small ruminant pregnancy; cattle and horse abdomen and thorax

LINEAR & MICROCONVEX

171 · Abdominal linear, 40 mm face — 6–14 MHz
Horse tendons; small-animal lungs and abdomen

150 / 195 · Abdominal microconvex, R15 mm — 4–9 MHz
Abdomen, AFAST and TFAST in dogs and cats; horse musculoskeletal; sheep and goat pregnancy

181 · Abdominal microconvex, R11 mm — 4–9 MHz
Dog and cat abdomen, AFAST/TFAST

OPU & BACKFAT

181/A · Microconvex for OPU, R11 mm — 4–9 MHz
OPU system (3.0 m cable)

174 · Backfat (pigs), 147 mm face — 2–5 MHz
Backfat and loin muscle depth in pigs

188/2 · 188/4 · Backfat (beef cattle), 180 mm face — 2–5 MHz
Backfat and loin muscle in beef cows (2.3 m or 4 m cable)

The set includes

- iScan 2 MULTI scanner
- External lithium-ion battery
- Battery charger
- USB storage adapter
- Carrying straps
- Transport case
- User guide

Optional accessories

- Wi-Fi antenna (images to mobile app)
- Stand
- OPU
- Extension and T-holder for rectal linear or rectal convex probe
- Shade cover / bag with shade cover
- Transport bag, suspenders, waist belt
- Standoff

Technical data

Imaging modes	B Mode, B+B Mode, B+M Mode
Frequency range	2–14 MHz (probe dependent)
Probes	150, 156, 161, 163, 171, 172, 174, 181, 181/A, 182, 188/2, 188/4, 191, 192, 195, 212, 225
Dynamic focus	Yes
Image management	Freeze, zoom 60–200%, full screen, invert left-right, image rotation, save images and cine loops
Imaging adjustment	Frame averaging, negative, LuciD
Post-processing	Yes (LuciD)
Presets	Save up to 5 presets
Measurements	Distance, follicle, area, area polygon, volume, grid, back-fat, needle; obstetric aging tables
Workflow tools	Personalized Quick Access Menu
Display	7″ LCD LED, 500 cd/m², 800 × 480 px
Controls	Membrane keyboard
Memory	200 images and 200 cine loops
Export	AVI, BMP · to USB flash drive · images to mobile app by Wi-Fi
Data port	6-pin connector for data transfer
Battery	Up to 7 h · charge up to 4 h
Start-up time	Under 20 s
Dimensions	25.6 × 16.5 × 6.9 cm
Weight	2.64 kg with probe and battery
Protection rating	IP32
Species	Cattle, horses, goats, sheep, pigs
Specialization	Reproduction, lungs

Technical service & support in Canada

Authorized service technicians

Our authorized service technicians maintain and repair your Draminski equipment here in Canada.

Loaner systems

If your scanner needs service, we can lend you a replacement system so your practice keeps working without interruption.

New trade-in program

Trade in your current ultrasound system and trade up to a newer Draminski scanner. Ask us for a trade-in evaluation.

Trade in. Trade up.

Ask us what your current system is worth toward an upgrade.

Flexible financing

Lease and financing options are available. Ask us for terms that fit your cash flow.

Manufacturer warranty

2 years on the scanner · 1 year on probes · 6 months on batteries

See scans & cine loops

Quotes, demos & service

ViiSona
ultrasound & education

SALES & TECHNICAL SERVICES
38 Rue Canvin
Kirkland, QC H9H 4S4
Canada

www.viisona.com · info@viisona.com · 1-800-672-9182

Starting prices in Canadian dollars (CAD), before applicable taxes and shipping. Final price depends on the probes and accessories selected. Specifications supplied by Dramiński S.A. and subject to change without notice. Specifications per draminski.com, October 2026. Pricing effective November 2026. Manufactured by Dramiński S.A., Poland · Sold and serviced in Canada by ViiSona.